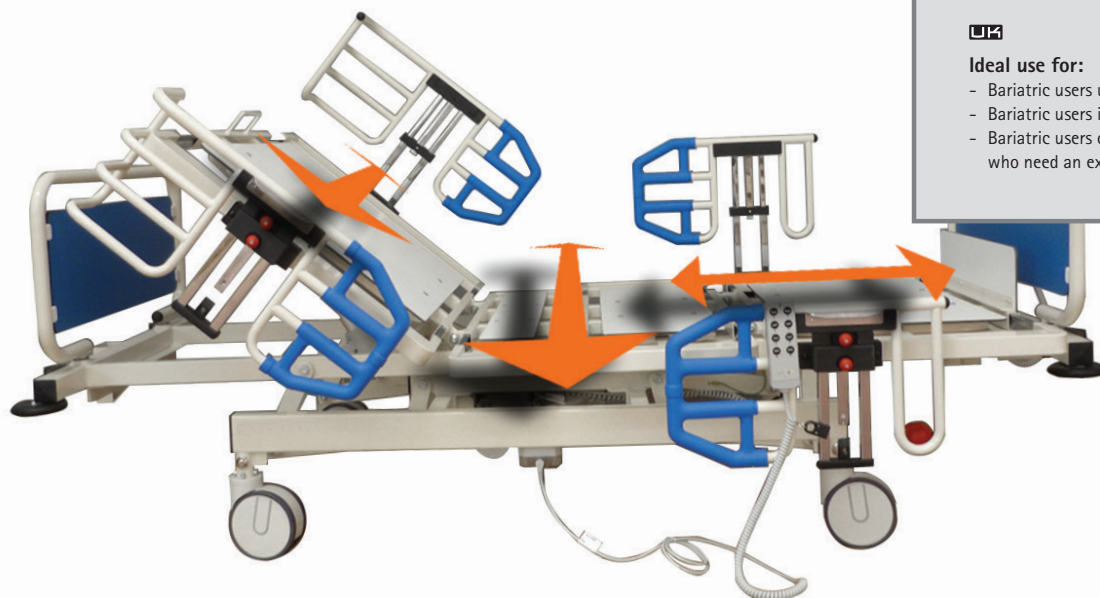
XXL-Rehab 400 Hospitalsseng

XXL-Rehab 400 bariatric hospital bed

DK

UK

MAX 400 KG

DK

Kan med fordel anvendes til:

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07/2014/JP

Instructions for using XXL-Rehab 400 bariatric hospital bed

⚠ CAUTION ⚠ CAUTION ⚠

- This product should not be used without proper instruction from a healthcare professional.
- Using the product for individuals weighing more than the maximum weight-bearing capacity may pose a risk for damaging the product.
- XXL-Rehab assumes no responsibility for any damage or injury caused by improper installation, assembly or use of this product.
- If components are damaged or missing, contact XXL-Rehab immediately. Do NOT use substitute parts.
- XXL-Rehab cannot be held liable for faults or accidents occurring after repairs by service staff without written authorisation from XXL-Rehab.

FOR USE

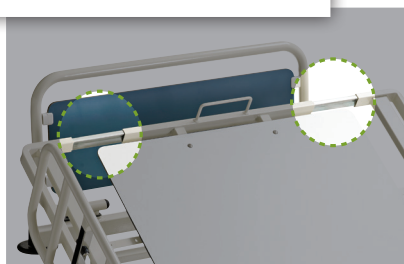
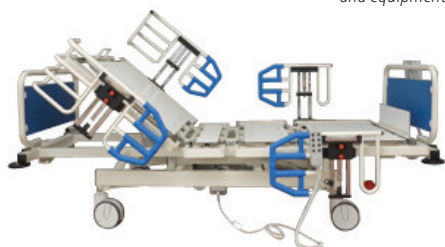
⚠ To ensure proper use of the XXL-Rehab 400 bariatric hospital bed, please read these instructions thoroughly before using the product. In doing so, pay particular attention to the warning notes marked with ⚠.

Also remember that in order to prevent the risk of injury to the user, the bed may be operated only by trained caregivers or aides or by individuals who have been instructed in its proper use.

INTENDED USE

The XXL-Rehab 400 bariatric hospital bed can be used for bariatric patients with a maximum weight of up to 400 kg. The bed is particularly suitable for hospitals and home care.

The bed can have a maximum load of 400 kg, including user and equipment.

Adjustable in width at both sides.

SPECIFICATIONS

Maximum weight-bearing capacity is	400 kg
Total length	231-251 cm
Total width	100-128 cm
Mattress size length	200 / 220 cm
Mattress size width	90 / 100 / 120 cm
Manually adjustable in length	200 / 220 cm
Manually adjustable in width	90 / 100 / 120 cm
Lying surface height	39-82 cm
Minimum clear distance from the floor	15 cm
Total weight	175 kg
Height adjustable	Yes
Separate control for caregiver e.g. Trendelenburg and Anti-Trendelenburg position	Yes
Electrical adjustment of backrest, leg section, thigh section, height of lying surface	Yes
Headboard up / down feature	Yes
Foot end up / down feature	Yes
Side rails removable without tools	Yes
Direction control on the wheels	Yes
Central brake	Yes
Magnetic lock on the user hand control	Yes
Built-in bed extension	Yes
Lock of individual electric functions performed from the caregiver hand control	Yes
Separate adjustment of backrest	Yes
Coupled adjustment of lying surface	Yes
Removable lifting pole	Yes
Removable foot gable	Yes
High slatsides	Yes
Frame	Powder-coated steel
Wheel diameter	145 mm
Ingress for water	66 IP
ISO - classification	ISO 18 12 10 - 30

ORDERING

Model	Order No
XXL-Rehab 400 bariatric hospital bed	2314-120-220
Accessories	
Side rails 2314-999-008	2314-999-008
Lifting pole 2314-999-010	2314-999-010
IV pole 2314-999-012	2314-999-012

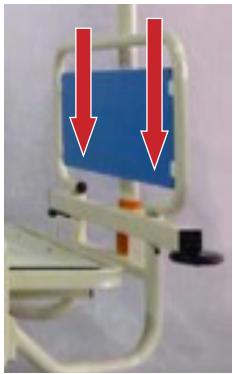
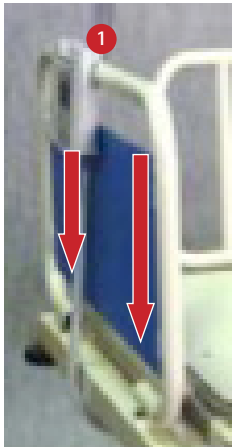
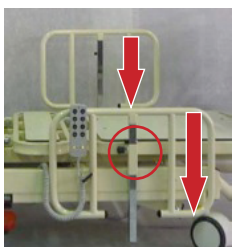
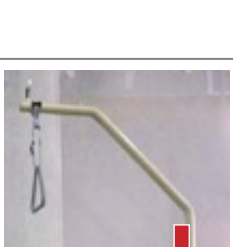
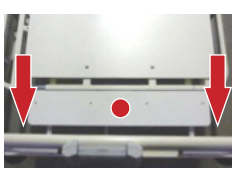
User hand control.

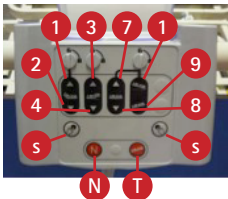
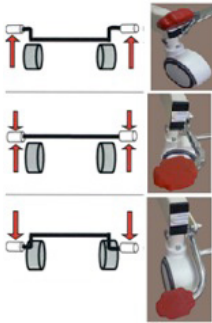
Caregiver control panel.

COMPONENT LIST TO ASSEMBLE

Assembling of head board	<p>Insert head board in the proper sockets on bed main frame, as red arrows show in the picture.</p> <p>TECHNICAL NOTES Head board sockets are designed to guarantee fast and reliable emergency removal.</p>	
Assembling of foot board	<p>Insert foot board in the proper sockets on bed main frame, as red arrows show in the picture.</p> <p>TECHNICAL NOTES Head board sockets are designed to guarantee fast and reliable emergency removal.</p>	
Caregiver control panel 1	<p>Attach the caregiver control panel to the upper part of the foot board.</p> <p>The caregiver control panel can quickly be removed in case of an emergency situation.</p>	
Adjustment of side rails	<p>Insert each side rail in its proper socket, as red arrows show in the picture. They can be adjusted at 3 different heights by a comfortable and reliable "click-clack" device as red circle shows in the picture: just pull the black knob, lift up the side rail, take it to the desired height, release the black knob and verify if the side rail is locked.</p>	
Assembling of lifting pole with IV hooks	<p>Insert bariatric lifting pole with IV hooks in its proper socket behind head board, on main frame as red arrow shows in the picture.</p> <p>TECHNICAL NOTES Bariatric lifting pole with IV hooks is made of 2 high-strength steel tubes.</p>	

Section widening – in width	<p>The bed can be adjusted in width from 90/100/110/120 cm. Sections can be adjusted to 2 different sizes by a comfortable and reliable "clickclack" device as shown in the picture: just pull the black knob with a hand and with the other hand pull out the section to the desired size, release the black knob and verify if the section is locked. Remember to adjust each section.</p>	
Bed extension – in length	<p>Turn the two black hand wheel knobs (right and left side) below the foot section and pull the foot section out, then extract the mattress holder from its socket and place it as shown in the picture.</p>	
Patient control panel	<p>Use of patient control panel can be disabled with the included magnetic key (M) by pulling it over "lock" symbol on the control panel (SZ).</p> <p>Lamp lights up: control panel is activated. Lamp off: control panel is disabled / locked.</p>	

COMPONENT LIST TO ASSEMBLE

Caregiver control panel S N T	<p>On the caregiver control panel there are 3 locking buttons above the electronic functions. By turning the locking knob in horizontal position the function of both patient and caregiver control panel is blocked.</p> <p>"S" button is activated by the caregiver to prevent that the electronic functions are used simultaneously.</p> <p>Red key "N": Emergency; brings the bed down.</p> <p>Red Key "T": Trendelenburg position; brings the bed in the above position.</p>	
Back rest (0°-54°) 1 2	<p>Push button 1 with the equivalent explanatory graphic to RAISE the BACKREST as shown in the picture (HANDSET). Push button 2 to LOWER the BACKREST.</p>	
Leg rest (0°-24°) 3 4	<p>Push button 3 with the equivalent explanatory graphic to RAISE the LEGREST as shown in the picture (HANDSET). Push button 4 to LOWER the LEGREST.</p>	
Back rest & leg rest simultaneously 5 6	<p>On Patient control panel push button 5 with the equivalent explanatory graphic to RAISE BACKREST & LEGREST SIMULTANEOUSLY as shown in the picture (HANDSET). On the patient control panel push button 6 to LOWER BACKREST & LEGREST SIMULTANEOUSLY.</p>	
Adjusting height (37-80 cm) 7 8	<p>Push button 7 with the equivalent explanatory graphic to RAISE THE BED as shown in the picture (HANDSET). Push button 8 to LOWER THE BED.</p>	
Trendelenburg position (max +/- 18°) 9	<p>Push button 9 with the equivalent explanatory graphic to get TRENDLENBURG position as shown in the picture (SUPERVISOR).</p>	
Continuous beep sound	<p>The operation is blocked or the function is completed.</p>	
Wheel steering	<p>A 3 step adjustable central braking pedal is placed on the steering wheels in the foot end of the bed.</p> <p>Steering position Push the pedal up until the steering castor is engaged. This adjustment allows driving the bed on a straight line.</p> <p>Swivelling position Push the pedal in the neutral position, castors are free to swivel, they are all unlocked.</p> <p>Braked position Push the pedal down, castors are all locked now.</p>	

MAINTENANCE (3 TIMES A YEAR) AND CLEANING

- After a week or two of usage, check all the bolts to be sure that none of them have worked loose and that all pins are in their normal location and securely fastened.
 - Check all welds.
 - Castors with brake: often check the brake blocks.
 - Often check that screws and nuts are kept well tightened down.
 - Use mild soap and water or a disinfectant to clean the bed.
- Wet a cloth with the cleaning product and wring dry.

⚠ Failure to inspect the bolts as directed may cause functional failures.

ACCESSORIES

Side rails, lifting pole and IV pole.

WARRANTY

The XXL-Rehab 400 bariatric hospital bed is covered by one year warranty from the invoice date for materials and workmanship.

XXL-Rehab guarantees that the product is flawless when it is delivered. If the product has a defect or is damaged when it is delivered, you should immediately contact XXL-Rehab.

Products which have been used for purposes other than the intended use; or which have been subject to negligence, abuse, improper storage or handling, improper operation, unauthorized modifications or damages beyond normal wear and tear as determined by XXL-Rehab, are not covered by this warranty.

If weight capacity is exceeded, the warranty will be void. Any unauthorized repairs to product / part, as well as tampering with any components, will void the warranty.

This warranty does NOT cover normal maintenance such as cleaning, adjustment, lubrication and updating of equipment or parts.

WARNINGS

Inform the patient of these rules of security:

- ⚠ The bed is only suitable for people over 12 years old.
- ⚠ The safe working load (mattress included) is 400 kg.
- ⚠ Persons who have not read the user manual must not use the bed.
- ⚠ It is not allowed to use the bed for other purposes than the ones stated in INTENDED USE.
- ⚠ For your safety, during bed transfers, be sure that the sections are in flat position and the height of the bed is at least 55 cm from the ground.
- ⚠ Remove any objects or obstacles that might be present in the working area around the bed.
- ⚠ Remember to lock the castors after moving the bed.
- ⚠ When not in use hook the hand control on the top of the side rails.
- ⚠ Make sure that all nuts, bolts and screws are kept well tightened.
- ⚠ For safety reasons, make sure to replace worn-out or damaged parts in due time.
- ⚠ Failure to inspect the bolts as directed may cause functional failures.

MANUFACTURER

XXL-Rehab
Fuglebækvej 1D
DK-2770 Kastrup
Denmark
Tel.: +45 70 25 25 22
Fax: +45 70 25 25 44
www.xxl-rehab.com

XXL-Rehab 400 bariatric hospital bed is CE marked in accordance with 93/42/CE directive (class 1) and CEI EN 60601-2-38; UNI EN 1970:2006; 2006/42/CE; CND: z12030799 – RDM: 106290/R – GMDN: 35563.

07/2014/JP

Strongly engineered and functionally designed exclusively for bariatric users!

XXL-Rehab Minimaxx

XXL-Rehab Minimaxx
Push Motor

XXL-Rehab Manual
Wheelchair - The Eclipse

XXL-Rehab Manual
Wheelchair - The Eclipse Tilt

XXL-Rehab
HD Electric Wheelchair

XXL-Rehab Wheelchair Cushions

XXL-Rehab Shower Bench

XXL-Rehab
Bedside commode

XXL-Rehab
Shower Commode

XXL-Rehab
Shower Commode Wheeler

XXL-Rehab
Shower Commode Rise N Recline

XXL-Rehab
Shower Commode Tilt

XXL-Rehab Barton Transferstol

XXL-Rehab Patient/Dinner chair

XXL-Rehab Rollator King

XXL-Rehab Rollator Support

XXL-Rehab Stand Tall

XXL-Rehab Bariatric walking frame

XXL-Rehab Canes

XXL-Rehab Multifunktional Scale

XXL-Rehab Lift and Slings

XXL-Rehab 400 bariatric hospital bed

XXL-Rehab Home Care Bed

XXL-Rehab Floor raiser

XXL-Rehab
Comfort Foam Mattress

XXL-Rehab
Low Air Loss mattress

XXL-Rehab
Alternating mattress

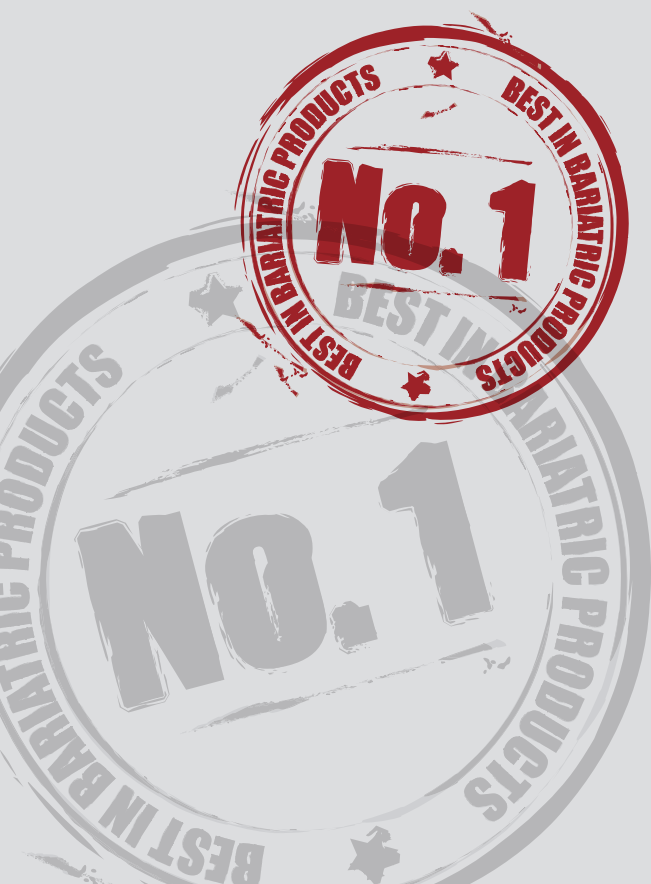
XXL-Rehab Bariatric suit

BRUGSANVISNING · USER MANUAL

mauryathome.dk

07-2014 XXL-Rehab 400 Hospitalsseng / XXL-Rehab 400 bariatric hospital bed

XXL REHAB®

XXL-Rehab · Fuglebækvej 1D · DK-2770 Kastrup · Denmark
Telephone: +45 7025 2522 · Fax: +45 7025 2544
info@xxl-rehab.com · www.xxl-rehab.com